



# Generation of a growth curve for iPS cells in a feeder-free culture by noninvasive image analysis.

Takayuki Uozumi<sup>1</sup>, Hiroaki Kii<sup>1</sup>, Momotaro Ishikawa<sup>1</sup>, Chieko Nakada<sup>1</sup>, Hirotada Watanabe<sup>1</sup>, Yasuyuki Asai<sup>2</sup>, Yasujiro Kiyota<sup>1</sup>  
<sup>1</sup> Instruments Company, NIKON CORPORATION, Yokohama, Japan, <sup>2</sup> ReproCELL, Inc, Yokohama, Japan

## Introduction

Maintaining the undifferentiated state of iPS cells is an essential, but difficult process to control. The process requires improvement in the following three areas:

1. Feeder-free culturing: To eliminate culture variability and to reduce culturing costs.
2. Automation: To reduce culturing costs, and to achieve consistency in the supply of iPS cells.
3. Establishment of a quantitative index: To quantify culturing status and provide a clear index for evaluation.

As a means for implementing these improvements, a system has been developed for time-lapse observation of feeder-free cultures within which iPS colonies can be automatically detected and measured.

The Nikon BioStation CT cell culture observation system was used to capture time-lapse images of the entire surface of a 6-well plate (100 images with a 2x objective (4x4mm)) over a 1 week period of feeder-free culturing (using ReproCELL ReproFF2), from seeding to passage. Once a culture dish is inserted and an observation schedule is set, the system automatically captures images while maintaining the observation position even with daily medium replacements, allowing the growth of each colony to be tracked with ease.

CL-Quant image analysis software was taught to detect iPS colonies, thus allowing iPS colonies to be identified within phase-contrast images. Changes within each colony were measured and plotted as a colony growth curve which could then be summarized per well, each representing a distinct culturing condition.

The system was used to generate growth curves for iPS cells cultured over a 1 week period (from seeding to passaging) while their undifferentiated state was maintained. Quantifying the rate of growth and comparing it to the undifferentiated state established a benchmark for this process.

## Methods


### iPS cell culture and growth



### iPS cell culture and growth

- iPS cells were generated from adult human dermal fibroblast (HDF) by retroviral transduction of four factors: Oct3/4, Sox2, Klf4, and c-Myc [1].
- Six days after transduction, the cells were harvested by trypsinization and plated onto feeder free cells using a 6-well plate. Repro FF2 medium ReproCELL was used.

### 1 week culture schedule

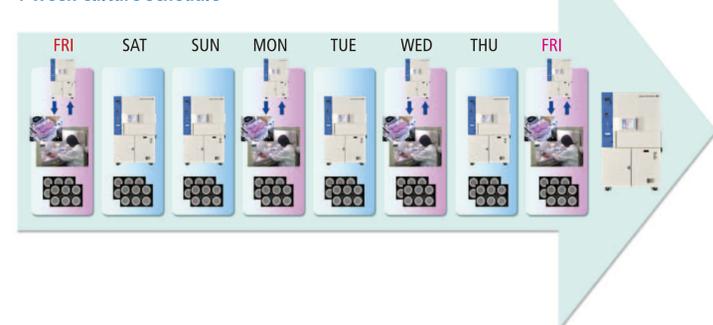

### Basic condition (1 week)

- Culture conditions  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6-well plate, Culture medium: ReproFF2, Medium change: Mon, Wed, Fri/ 1 week
- Capture conditions  
Magnification: x2 (Full scan), Culture day: 1 week, Imaging interval: 3 or 12 hours

### Image acquisition


- Phase contrast images of a 6-well plate were acquired using BioStation CT (NIKON CORPORATION) at 37°C, 5% CO<sub>2</sub>.
- High-precision repeatability: Accurate tracing of same cells, as well as X-Y positions (±5µm) for each vessel.
- Full Scan mode (2x) was set for acquiring full-well tiling images.

### Image analysis

- Full-well tiling images were loaded into CL-Quant software (NIKON CORPORATION) and stitched to create 10,000-by-10,000 pixel composite images for analysis.
- All image analyses (colony segmentation, measurement etc.) were performed using functions built into CL-Quant software [2].

## Algorithm development


### Colony detection & calculation flow

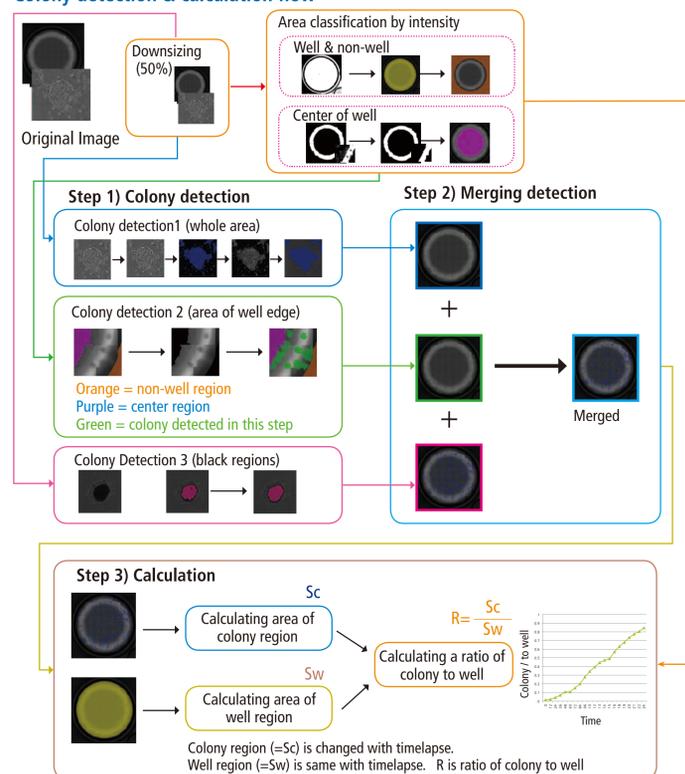

## Confirmation of algorithm

### Algorithm development

Nikon has developed an algorithm that enables temporal quantitative analysis of colony proliferation, from immediately after cell seeding to cell confluency.

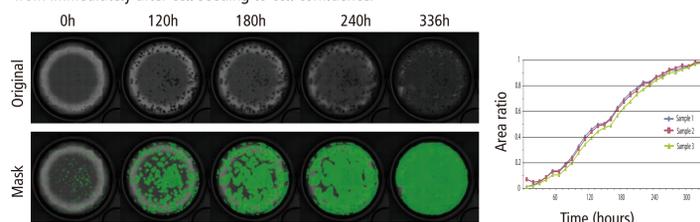


### Experiment procedures

- Culture conditions  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6 well plate, Culture medium: ReproFF2  
Medium change: Mon, Wed, Fri/ 1 week
- Capture conditions  
Magnification: x2 (Full scan), Culture day: 2 weeks, Imaging interval: 12hours

### Comparison of results from detection using algorithm and manual detection

Manual versus auto-detection with CL-Quant calculations were made. The results were similar from center to edge of dish.

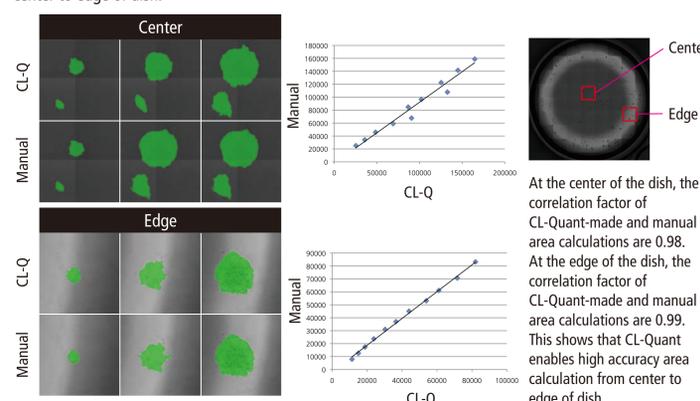


### Experimental procedure

- Culture condition  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6 well plate, Culture medium: ReproFF2  
Medium change: Mon, Wed, Fri/ 1 week
- Capture condition  
Magnification: x2 (Full scan), Culture day: 1 week, Imaging interval: 3hours

## Results

### Experiment 1) Using the CL-Quant algorithm for standard cell culture

CL-Quant enables image analysis of standard 1 week iPS cell-culture sequence.

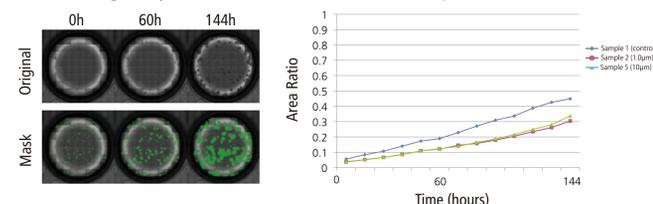


### Experiment procedures

- Culture conditions  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6 well plate, Culture medium: ReproFF2  
Medium change: Mon, Wed, Fri/ 1 week
- Capture conditions  
Magnification: x2 (Full scan), Culture day: 1 week, Imaging interval: 12hours

### Experiment 2) Comparison of cell seeding density

The results indicate that growth of colonies depends on seeding density (for 1 week cell cultivation).

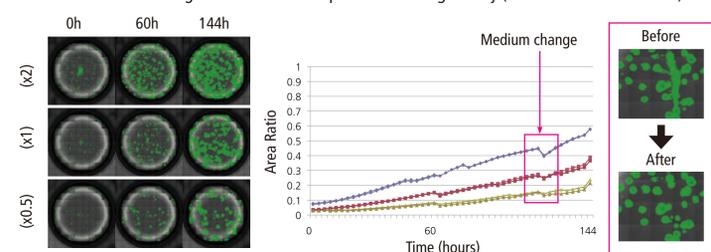


### Experiment procedures

- Culture conditions  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6 well plate, Culture medium: ReproFF2  
Medium change: Mon, Wed, Fri/ 1 week
- Capture conditions  
Magnification: x2 (Full scan), Culture day: 2 weeks, Imaging interval: 3 hours

### Experiment 3) Drug addition

CL-Quant quantifies the influence of Tzv-addition, which affects a cell's survival rate.

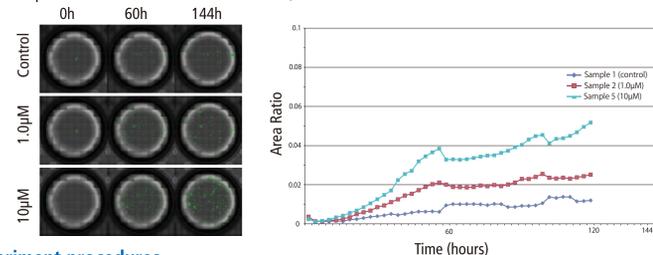

### Experiment procedures

- Culture condition  
Cell: hiPSs:201B7, Coating: Matrigel, Culture vessel: 6 well plate, Culture medium: ReproFF2  
Medium change: Mon, Wed, Fri/ 1 week
- Capture conditions  
Magnification: x2 (Full scan), Culture day: 5 days, imaging interval: 3 hours

### Thiazovinin Mechanism Overview

#### Description

Tzv is a novel ROCK inhibitor and Rho-ROCK axis regulates cell-ECM and cell-cell adhesion. Thiazovinin protects human embryonic stem cells (hESCs) in the absence of ECM by regulating E-cadherin mediated cell-cell interaction [3]. This observation suggests that Thiazovinin promotes cell survival. Thiazovinin, in combination with inhibitors of the TGF-β receptor and MEK pathway, has shown to improve reprogramming efficiency by more than 200-fold [4].

## Reference

1. Takahashi, K., Yamanaka, S. Induction of pluripotent stem cells from adult human fibroblasts by defined factors. Cell. (2007) Nov 30;131(5):861-72.
2. Alworth, SV, Watanabe, H. & Lee, JSJ. Teachable, High-Content Analytics for Live-Cell, Phase Contrast Movies. J Biomol Screening. (2010) 8, 968-977.
3. Xu, Y., Ding, S. Revealing a core signaling regulatory mechanism for pluripotent stem cell survival and self-renewal by small molecules. Proc Natl Acad Sci USA. (2010) 107: 8129-8134.
4. Lin, T., Ding, S. A chemical platform for improved induction of human iPS cells. Nature Methods (2009) 6: 805-808. (URL:https://www.stemgent.com/products/show/72)

## Conclusion

1. Nikon has developed an algorithm that enables quantitative analysis of colony proliferation, from immediately after cell seeding to confluency. Comparing manual and algorithm-detection, the correlation factors in full image areas have similar results: 0.98 at the center of the dish, and 0.99 at edge.
  2. CL-Quant enables analysis of effects under various conditions.
    - Colony area proliferation depends on cell seeding density.
    - CL-Quant quantifies the effects of thiazovinin density.
- CL-Quant has achieved non-invasive quantification of change of culture conditions of iPS cells in a feeder-free culture, based on the ratio of colony/bottom of dish area.